
drawings 1-7

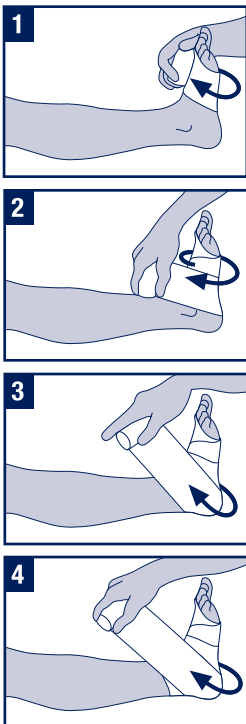Elodur® forte:

94% Baumwolle/cotton
6% Elasthan/elastane

Eloflex®:

57% Baumwolle/cotton
33% Viskose/viscose
8% Poliamid/polyamide
2% Elasthan/elastane

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Elodur® forte

EN

Medium to long stretch bandages for support, relief and compression.

Dressing instruction:

We recommend that bandaging be performed by a trained person.

Indications:

For support and relief in diseases of the ligament and postural system and for sports injuries, e.g. distortions, contusions, luxations, tendovaginitis, rhexes of the muscle fibre and post-traumatic oedemas/haematomas. For compression dressing for phlebological indications, e.g. the follow-up treatment of venous diseases.

Eloflex®: prophylaxis of thrombosis.

Contraindications:

Advanced peripheral arteriosclerosis, decompensated cardiac insufficiency, septic phlebitis, phlegmasia coerulea dolens, incompatibility to fabric.

Cautions:

- Weeping dermatosis
- Cutaneous infections
- Open wounds must be covered with an appropriate dressing before compression bandaging is put on.
- Severely compromised skin sensibility and impaired sensitivity of the limb
- Advanced peripheral neuropathy
- Conditions in which increased venous and lymphatic return is not desired.

The attending physician decides whether compression in general can be applied to patients suffering from the following diseases: rheumatoid arthritis, complex regional pain syndrome, CRPS (M. Sudeck), malignant lymphedema, gangrene.

Longitudinal elasticity:

Elodur® forte: ca. 175%
Eloflex®: ca. 180%

Special notes:

Due to their high resting pressure, long-stretch bandages should be removed at night (does not apply to **Eloflex®**-prophylaxis of thrombosis).

Any serious incident that has occurred in relation to this device should be reported to BSN medical Inc. (BSN medical GmbH within the European Union) and the competent authority of your state.

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