

# Cutimed® Siltec® Plus Heel

Cutimed® Siltec® Plus Heel is a silicone coated, sterile, absorbent polyurethane foam dressing for atraumatic dressing changes that contains super absorbent stripes which absorb and lock wound exudate.

## Product indication

The Cutimed® Siltec® Plus Heel assortment is indicated for exuding wounds such as: venous and arterial ulcers, pressure ulcers, diabetic foot ulcers surgical incisions, skin grafts and donor sites, lacerations or abrasions. It is recommended for wounds with low to high exudate levels.

## Product benefits



### Reliable exudate management

The vertical absorption makes sure exudate is taken up reliable and kept inside the dressing - even under compression. Super absorbent stripes additionally support the absorption capacity.<sup>1,4</sup>



### Gentle with high adherence

The perforated silicone wound contact layer - with higher adherence than Cutimed® Siltec® Heel 3D - prevents adherence to newly formed tissue, protects wound margins and ensures atraumatic, pain-reduced dressing changes.<sup>1,2,3</sup>



### Water-repellent and breathable backing

The outer film is bacteria proof, water-repellent and ensures proper fluid evaporation for changing stages of saturation.



### Cushioning effect

This dressing reduces pressure peaks. Its cushioning properties provide good wearing comfort.




### Special shape

Especially designed to fit the heel and provide good conformability. Also adaptable to similar body parts, such as elbows or knees. The dressing can also be cut to size.<sup>1,2</sup>



## Applications

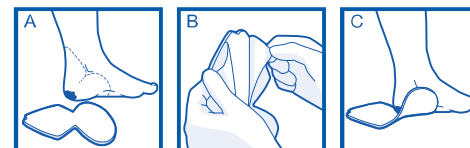

1 Seckam A., A multicentre, observational evaluation of the product characteristics of two absorbent foam dressings. *Br J Nurs.* Jun 27, 2019;28(12):S10-S17.

2 King B., A clinical evaluation of 20 patients when using a new absorbent silicone foam wound dressing: Cutimed Siltec B. *Wounds UK.* 2018; Vol 14. No 3.

3 Zimmermann S., *In vitro* exudate management performance of a new hydrophilic foam dressing range. *Poster EWMA.* 2017; Amsterdam.

4 Bateman SD., 150 patient experiences with a soft silicone foam dressing. *Br J Nurs.* Jun 25-Jul 8, 2015;24(12):S16, S18-23.

# Cutimed® Siltec® Plus Heel

Ordering Information	Size	Product Code	Pack size	Nappi Code
----------------------	------	--------------	-----------	------------

## Cutimed® Siltec® Plus Heel

Highly absorbent  
Silicone foam dressing  
with stronger initial  
adherence for the heel

16 x 24 cm

73288-05

5

274590-001



LATEX  
**FREE**  
FORMULA

## Icon overview

### Improved guidance with symbols

Exudate	Dry Low Moderate High Excessive Low to Excessive
Wound depth	superficial deep superficial and deep
Wound state	Infected Sloughy Granulating Epithelializing Necrotic
Both side application	Not on pack (in IFU only)
Secondary dressing needed	

Find icons here



# Cutimed® Siltec® Heel 3D

Cutimed® Siltec® Heel 3D is a silicone coated, sterile, absorbent polyurethane foam dressing for atraumatic dressing changes that contains super absorbent stripes which absorb and lock wound exudate.

## Product indication

Cutimed® Siltec® Heel 3D is indicated for exuding wounds such as: venous and arterial ulcers, pressure ulcers, diabetic foot ulcers, surgical incisions, skin grafts and donor sites, lacerations or abrasions.

## Product benefits



### Reliable exudate management

The vertical absorption makes sure exudate is taken up reliably and kept inside the dressing - even under compression. Super absorbent stripes additionally support the absorption capacity.<sup>1,4</sup>



### Gentle to wound bed

The skin-friendly, perforated silicone wound contact layer prevents adherence to newly formed tissue, protects wound margins and ensures atraumatic, pain-reduced dressing changes.<sup>1,2,3</sup>



### Water-repellent and breathable backing

The outer film is bacteria proof, water-repellent and ensures proper fluid evaporation for changing stages of saturation.



### Cushioning effect

This dressing reduces pressure peaks. Its cushioning properties provide good wearing comfort.




### Special shape

Especially designed to fit the heel and provide good conformability. The dressing can also be cut to size.<sup>1,2</sup>



## Applications

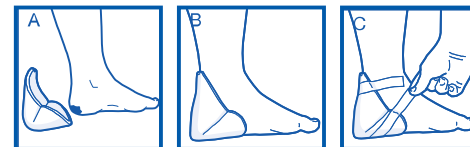

1 Seckam A., A multicentre, observational evaluation of the product characteristics of two absorbent foam dressings. *Br J Nurs.* Jun 27, 2019;28(12):S10-S17.

2 King B., A clinical evaluation of 20 patients when using a new absorbent silicone foam wound dressing: Cutimed Siltec B. *Wounds UK.* 2018; Vol 14. No 3.

3 Zimmermann S., *In vitro* exudate management performance of a new hydrophilic foam dressing range. *Poster EWMA.* 2017; Amsterdam.

4 Bateman SD., 150 patient experiences with a soft silicone foam dressing. *Br J Nurs.* Jun 25-Jul 8, 2015;24(12):S16, S18-23.

# Cutimed® Siltec® Heel 3D

Ordering Information	Size	Product Code	Pack size	Nappi Code
----------------------	------	--------------	-----------	------------

**Cutimed® Siltec®  
Heel 3D**

Highly absorbent  
Silicone foam dressing  
for the heel, gently  
adherent

16 x 24 cm

73286-01

5

274591-001



**LATEX  
FREE  
FORMULA**



## Icon overview

### Improved guidance with symbols

Exudate	Dry Low Moderate High Excessive Low to Excessive
Wound depth	superficial deep superficial and deep
Wound state	Infected Sloughy Granulating Epithelializing Necrotic
Both side application	Not on pack (in IFU only)
Secondary dressing needed	

Find icons here

