

Cutimed® Siltec® Sacrum

Cutimed® Siltec® Sacrum is a silicone coated, sterile, absorbent polyurethane foam dressing for atraumatic dressing changes that contains super absorbent stripes which absorb and lock wound exudate. It has an additional silicone adhesive border, that allows secure and gentle fixation of the dressing.

Product indication

The Cutimed® Siltec® Sacrum assortment is indicated for exuding wounds such as: venous and arterial ulcers, pressure ulcers, diabetic foot ulcers, surgical incisions, skin grafts and donor sites, lacerations or abrasions. Cutimed® Siltec® Sacrum is recommended for wounds with low to high exudate levels.

Product benefits



Reliable exudate management

The vertical absorption makes sure exudate is taken up reliable and kept inside the dressing - even under compression. Super absorbent stripes additionally support the absorption capacity.¹



Gentle to wound bed

The skin-friendly, perforated silicone wound contact layer prevents adherence to newly formed tissue, protects wound margins and ensures atraumatic, pain-reduced dressing changes.³



Self-adhesive silicone border

The self-adhesive silicone border, with low allergy potential, keeps the dressing in place and allows an easy application.³



Shower-proof

The highly breathable, water-repellent backing protects the dressing and even makes it possible to take a shower without removing or changing the dressing.




Special shape

Especially designed to fit the sacral area and provide good conformability. The dressing can also be cut to size.^{1,2}



Applications

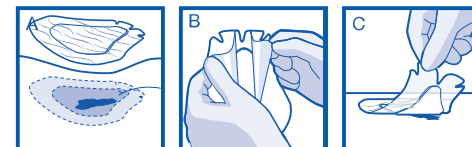

¹ Seckam A., A multicentre, observational evaluation of the product characteristics of two absorbent foam dressings. *Br J Nurs.* Jun 27, 2019;28(12):S10-S17.

² King B., A clinical evaluation of 20 patients when using a new absorbent silicone foam wound dressing: Cutimed Siltec B. *Wounds UK.* 2018; Vol 14, No 3.


³ Zimmermann S., *In vitro* exudate management performance of a new hydrophilic foam dressing range. *Poster EWMA.* 2017; Amsterdam.

Cutimed® Siltec® Sacrum

Ordering Information	Size	Product Code	Pack size	Nappi Code
Cutimed® Siltec® Sacrum Highly absorbent Silicone foam dressing with silicone adhesive edge for the sacral area	17,5 x 17,5 cm	73287-00	5	274588-001
	23 x 23 cm	73287-01	5	274589-001



LATEX FREE FORMULA

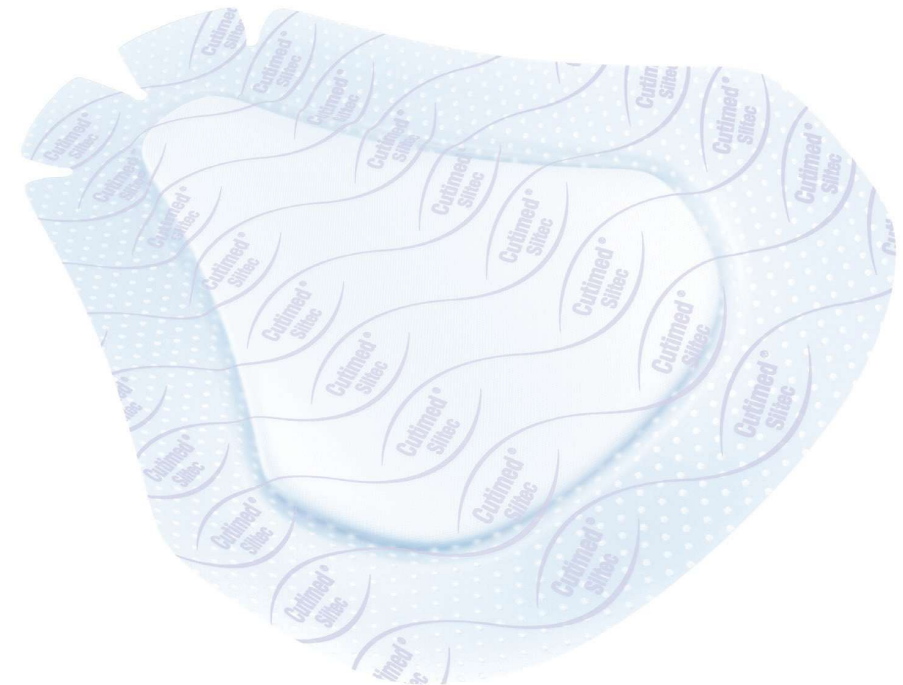

Icon overview

Improved guidance with symbols	
Exudate	Dry Low Moderate High Excessive Low to Excessive
Wound depth	superficial deep superficial and deep
Wound state	Infected Sloughy Granulating Epithelializing Necrotic
Both side application	Not on pack (in IFU only)
Secondary dressing needed	

Find icons here

